For Research Use Only

CoraLite® Plus 488-conjugated Cytokeratin 19 Polyclonal antibody



Catalog Number: CL488-14965 Featured Product

Basic Information

Catalog Number: CL488-14965 Concentration:

1000 ug/ml Source:

Rabbit Isotype:

Immunogen Catalog Number:

AG6830

GenBank Accession Number:

BC002539 GeneID (NCBI): 3880

UNIPROT ID: P08727

Full Name: keratin 19 Calculated MW:

44 kDa

Purification Method:

Antigen affinity purification

Recommended Dilutions: IF/ICC: 1:50-1:500

FC (Intra): 0.40 ug per 10⁶ cells in a

100 µl suspension

Excitation/Emission maxima

wavelengths: 493 nm / 522 nm

Applications

Tested Applications: IF/ICC, FC (Intra) Species Specificity: human, mouse, rat

Positive Controls:

IF/ICC: HeLa cells, FC (Intra): MCF-7 cells,

Background Information

Keratins are a large family of proteins that form the intermediate filament cytoskeleton of epithelial cells, which are classified into two major sequence types. Type I keratins are a group of acidic intermediate filament proteins, including K9-K23, and the hair keratins Ha1-Ha8. Type II keratins are the basic or neutral courterparts to the acidic type I keratins, including K1-K8, and the hair keratins, Hb1-Hb6. KRT19, one of type I keratins, is specifically expressed in the periderm, the transiently superficial layer that envelopes the developing epidermis. Due to its high sensitivity, KRT19 is the most used marker for detection of tumor cells disseminated in lymph nodes, peripheral blood, and bone marrow of breast cancer patients. This antibody, generated against full length KRT19 protein, also recognizes KRT17, another type I keratin homologous to KRT19.

Storage


Store at -20°C. Avoid exposure to light. Stable for one year after shipment.

Storage Buffer


PBS with 50% glycerol, 0.05% Proclin300, 0.5% BSA, pH7.3

Aliquoting is unnecessary for -20°C storage

Selected Validation Data

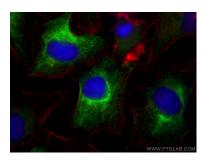

Immunofluorescent analysis of (-20°C Ethanol) fixed HeLa cells using CoraLite® Plus 488 Cytokeratin 19 antibody (CL488-14965) at dilution of 1:200, CL594-phalloidin (red).

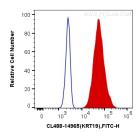

1X10^6 MCF-7 cells were intracellularly stained with 0.4 ug CoraLite® Plus 488 Anti-Human Cytokeratin 19 (CL488-14965) (red), or 0.4 ug Control Antibody. Cells were fixed with 4% PFA and permeabilized with Flow Cytometry Perm Buffer (PF00011-C).