

For Research Use Only

# CoraLite®594-conjugated p115, USO1 Polyclonal antibody



Catalog Number:CL594-13509

Featured Product

## Basic Information

Catalog Number:

CL594-13509

Size:

1000 µg/ml

Source:

Rabbit

Isotype:

IgG

Immunogen Catalog Number:

AG4431

GenBank Accession Number:

BC032654

GeneID (NCBI):

8615

UNIPROT ID:

O60763

Full Name:

USO1 homolog, vesicle docking protein (yeast)

Calculated MW:

962 aa, 108 kDa

Purification Method:

Antigen affinity purification

Recommended Dilutions:

IF/ICC 1:50-1:500

Excitation/Emission maxima wavelengths:

588 nm / 604 nm

## Applications

Tested Applications:

IF/ICC

Species Specificity:

human

Positive Controls:

IF/ICC : HeLa cells,

## Background Information

p115, also known as USO1, TAP (transcytosis-associated protein) or VDP (vesicle docking protein) is a general vesicular transport factor and plays an important role at different steps of vesicular transport. It is a 962-residue peripheral membrane protein which recycles between the cytosol and the Golgi apparatus during interphase (PMID: 9478999). p115 forms stable homodimers (PMID: 19247479). Rab1 recruits p115 to coat protein complex II (COPII) vesicles during budding from the endoplasmic reticulum, where p115 interacts directly with a select set of SNARE proteins (PMID: 10903204). p115 is required for intra-Golgi transport, and also functions in endoplasmic reticulum to Golgi trafficking, Golgi biogenesis and exocytotic transport (PMID: 19247479).

## Storage

Storage:

Store at -20°C. Avoid exposure to light.

Storage Buffer:

PBS with 50% Glycerol, 0.05% Proclin300, 0.5% BSA, pH 7.3.

Aliquoting is unnecessary for -20°C storage

For technical support and original validation data for this product please contact:

T: 4006900926


E: [Proteintech-CN@ptglab.com](mailto:Proteintech-CN@ptglab.com)

W: [ptgcn.com](http://ptgcn.com)

This product is exclusively available under Proteintech Group brand and is not available to purchase from any other manufacturer.

---

## Selected Validation Data

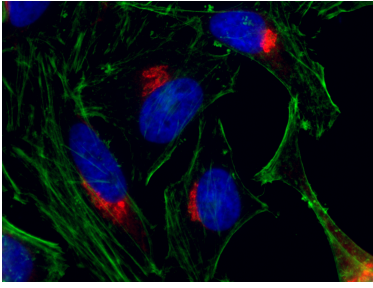

Immunofluorescent analysis of (4% PFA) fixed HeLa cells using CoraLite®594 p115, USO1 antibody (CL594-13509) at dilution of 1:200, 488-phalloidin (green).